






for Fast Blocking in 5 min

Bullet Blocking One for Western Blotting

- Features
- Fast blocking in **5 minutes.**
 - Ready to use.

Comparison Data

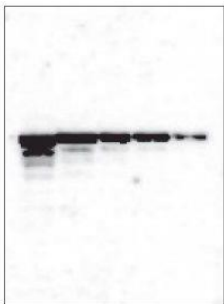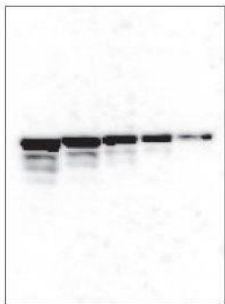Blocking Reagents	Bullet Blocking One	Blocking One	5% Skim Milk
Blocking Time	5 min.	30 min.	60 min.
Blotting Image			

* Blocking One is our standard blocking reagent.

(Conditions)

Blotting: Semi-dry Blotting Buffer Solution (Cat. No.: 30650) at 10 V for 30 min
 Washing: 0.1% t-TBS (Cat. No.: 12750)
 Blocking: Refer to the figure above
 1st antibody: Anti-Vimentin (C-20) (Rabbit) (Cat. No.: SC-7557-R) diluted 1:2,000, 1 hr. at RT
 2nd antibody: Anti-Rabbit IgG-HRP (Cat. No.: SC-2004) diluted 1:100,000, 1 hr. at RT
 Detection: Chemi-Lumi One (Cat. No.: 11644), 5 min. reaction time
 Detector: LAS-3000 (High mode), 15 min. exposure time
 Sample: 20 µg of HeLa cell extract, 5 serial two-fold dilution series

Ordering Information

Product Name	PKG size	Cat. No.	Price	Special Offer
Bullet Blocking One for Western Blotting	50 mL	13779-56	\$25.00	\$12.50
	200 mL (50 mL x 4)	13779-14	\$60.00	\$30.00
	1 L	13779-01	\$180.00	

Promo Code: **BB2022** Valid until **08/31/2022**






Compatible with Bio-Rad Mini-PROTEAN CosmoPAGE Precast Gel

Features


- Consistent performance
- Compatible with Bio-Rad Mini-PROTEAN Core
- Sharp, straight bands
- Long shelf life; Up to 12 months

Consistent Performance

CosmoPAGE gels perform high protein resolution and reproducibility even stored at ambient temperature for 50 days.

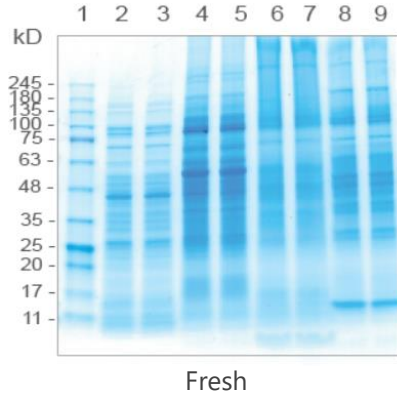

Fresh

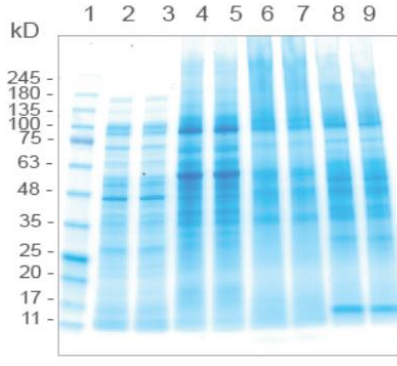

50 days at 37°C

(Conditions)
Gel : CosmoPAGE TG Gel 4-15 %
Voltage: 150V for an hour
Sample:
Lane 1. Prestained Protein Ladder, Broad Range, 3 µl (#NU01002)
Lane 2. 3. E. coli lysate
Lane 4. 5. Hek293T lysate
Lane 6. 7. Mouse brain membrane
Lane 8. 9. Mouse brain cytosol (10 µg each)

Ordering Information

Product Name	Cat. No.	PKG size	Price	Special Offer
CosmoPAGE TG (Tris-Glycine) Precast Gel 10%, 15 wells	NU001015	10 gels	\$89.00	\$62.30
CosmoPAGE TG (Tris-Glycine) Precast Gel 4-15%, 15 wells	NU041515	10 gels	\$89.00	\$62.30
CosmoPAGE Bis-Tris Precast Gel 12%, 15 wells	NU001215	10 gels	\$89.00	\$62.30
CosmoPAGE Bis-Tris Precast Gel 4-12%, 15 wells	NU041215	10 gels	\$89.00	\$62.30
CosmoPAGE MOPS Running Buffer, Powder	NU01004	5 x 1 L	\$25.00	



Promo Code: CP2022 Valid until 08/31/2022