






for Fast Blocking in 5 min

# Bullet Blocking One for Western Blotting

- Features
- Fast blocking in **5 minutes.**
  - Ready to use.

Comparison Data

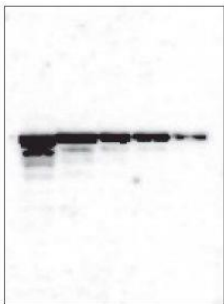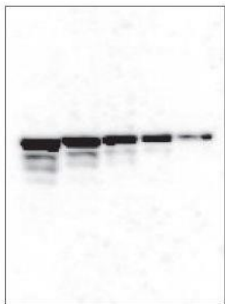Blocking Reagents	Bullet Blocking One	Blocking One	5% Skim Milk
Blocking Time	5 min.	30 min.	60 min.
Blotting Image			

\* Blocking One is our standard blocking reagent.

(Conditions)

Blotting: Semi-dry Blotting Buffer Solution (Cat. No.: 30650) at 10 V for 30 min  
 Washing: 0.1% t-TBS (Cat. No.: 12750)  
 Blocking: Refer to the figure above  
 1st antibody: Anti-Vimentin (C-20) (Rabbit) (Cat. No.: SC-7557-R) diluted 1:2,000, 1 hr. at RT  
 2nd antibody: Anti-Rabbit IgG-HRP (Cat. No.: SC-2004) diluted 1:100,000, 1 hr. at RT  
 Detection: Chemi-Lumi One (Cat. No.: 11644), 5 min. reaction time  
 Detector: LAS-3000 (High mode), 15 min. exposure time  
 Sample: 20 µg of HeLa cell extract, 5 serial two-fold dilution series

Ordering Information

Product Name	PKG size	Cat. No.	Price	Special Offer
Bullet Blocking One for Western Blotting	50 mL	13779-56	\$25.00	<b>\$12.50</b>
	200 mL (50 mL x 4)	13779-14	\$60.00	<b>\$30.00</b>
	1 L	13779-01	\$180.00	

Promo Code: BB22

Valid until 02/28/2022




# Prestained Protein Ladder



Cat. No.: NU01002

Prestained Protein Ladder, Broad Range (10-245 kDa)

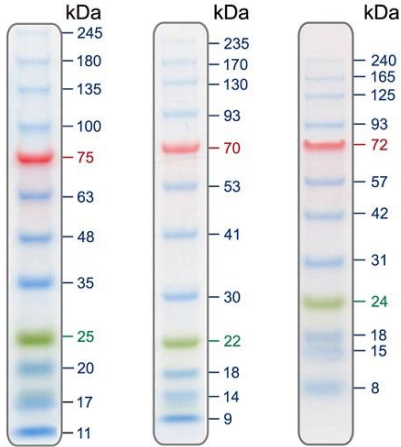


4~15% TG  
SDS buffer

4~12% Bis-Tris  
MOPS buffer

4~12% Bis-Tris  
MES buffer

Cat. No.: NU01009

Prestained Protein Ladder, Mid Range (15-180 kDa)

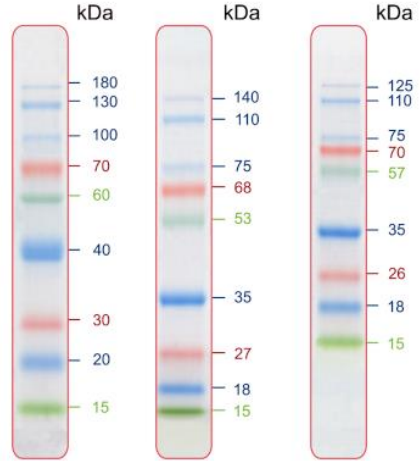

4~20% TG  
SDS buffer

4~12% Bis-Tris  
MOPS


4~12% Bis-Tris  
MES buffer

# Unstained Protein Ladder



Cat. No.: NU01010

Unstained Protein Ladder, Broad Range (10-200 kDa)

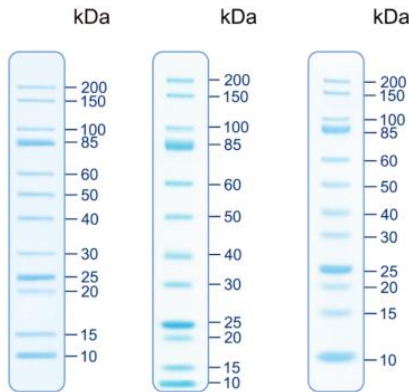

4~20% TG  
SDS buffer

4~12% Bis-Tris  
MOPS

4~12% Bis-Tris  
MES buffer

Ordering  
Information

Product Name	PKG size	Cat. No.	Price	Special Offer
Prestained Protein Ladder, Broad Range	500 µl	NU01002	\$ 100.00	<b>\$60.00</b>
Prestained Protein Ladder, Mid Range	500 µl	NU01009	\$ 100.00	<b>\$60.00</b>
Unstained Protein Ladder, Broad Range	500 µl	NU01010	\$ 100.00	<b>\$60.00</b>

Promo Code: BB22P

Valid until 02/28/2022

