

**20% OFF SALE**  
Valid till Aug 31, 2021

# CHEMILUMINESCENT HRP SUBSTRATES CHEMI-LUMI ONE




Get our Blocking Buffer for FREE!

For more detail please see back of the page.



## Chemi-Lumi One Series for HRP

Chemi-Lumi One is a series of high sensitive luminol-based chemiluminescence assay kits for Western Blotting. We offer three types of chemiluminescent substrates for Western blot detection with horseradish peroxidase enzyme (HRP).

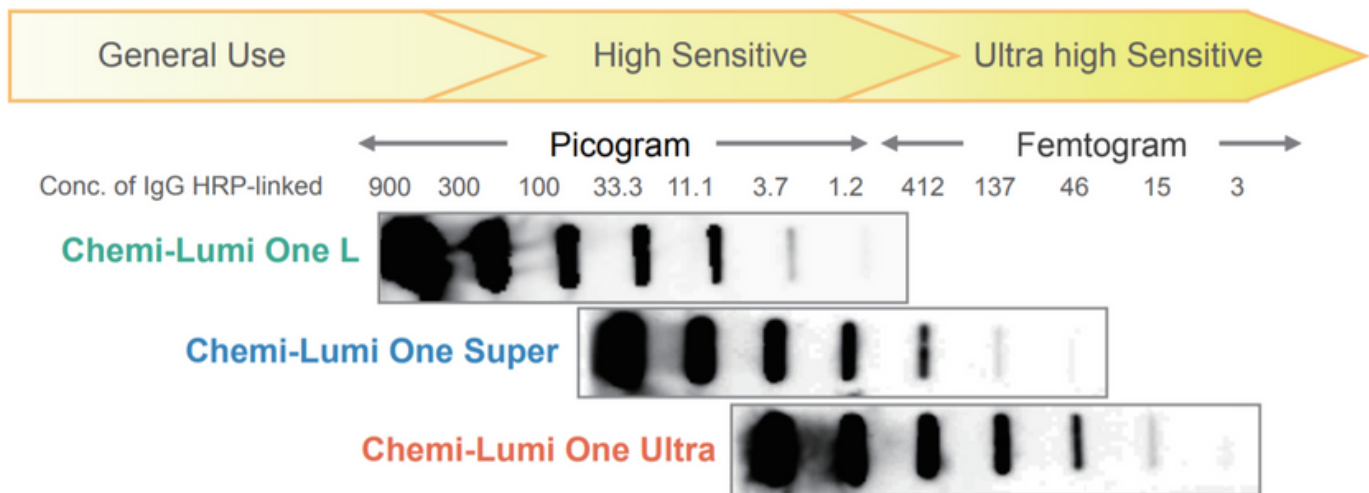

**Contact:**

TEL: 858-404-0403  
Email: info@nacalaiusa.com  
www.nacalaiusa.com

**Nacalai USA**  
Innovations for Life Sciences.



## PROMOTION CODE

**ChemiL20** 20%-OFF Chemi-Lumi One L, Super, Ultra




**ChemiL20-FB** 20%-OFF Chemi-Lumi One L, Super, Ultra + FREE 50 mL of Bullet Blocking Buffer

### Price List

Product Name		Storage	Cat. No.	Quantity	List Price	Price (US\$)
Chemi-Lumi One L	Luminol 250 mL, Peroxidase 250 mL Sufficient For: 4,000 cm <sup>2</sup> of membrane	4°C	07880-70	1 Kit	\$190.00	<b>\$152.00</b>
Chemi-Lumi One Super	Luminol 50 mL, Peroxidase 50 mL Sufficient For: 1,000 cm <sup>2</sup> of membrane	4°C	02230-30	1 Kit	\$200.00	<b>\$160.00</b>
	Luminol 10 mL, Peroxidase 10 mL Sufficient For: 200 cm <sup>2</sup> of membrane	4°C	02230-14	1 Kit (Trial Kit)	\$80.00	<b>\$64.00</b>
Chemi-Lumi One Ultra	Luminol 50 mL, Peroxidase 50 mL Sufficient For: 1,000 cm <sup>2</sup> of membrane	4°C	11644-40	1 Kit	\$298.00	<b>\$238.40</b>
	Luminol 10 mL, Peroxidase 10 mL Sufficient For: 200 cm <sup>2</sup> of membrane	4°C	11644-24	1 Kit (Trial Kit)	\$110.00	<b>\$88.00</b>

## Bullet Blocking One for Western Blotting

- Fast blocking in 5 minutes.
- Ready to use.

Blocking Reagents	Bullet Blocking One	Blocking One	5% Skim Milk
Blocking Time	5 min.	30 min.	60 min.
Blotting Image			

Sample: 20 µg of HeLa cell extract, 5 serial two-fold dilution series

Blotting: Semi-dry Blotting Buffer Solution (Product No. 30650) at 10 V for 30 min

Washing: 0.1% t-TBS (Product No. 12750)

Blocking: Refer to the figure above

1st antibody: Anti-Vimentin (C-20) (Rabbit) (Product No. SC-7557-R) diluted 1:2,000, 1 hr. at RT

2nd antibody: Anti-Rabbit IgG-HRP (Product No. SC-2004) diluted 1:100,000, 1 hr. at RT

Detection: Chemi-Lumi One (Product No. 11644), 5 min. reaction time

Detector: LAS-3000 (High mode), 15 min. exposure time

### Price List

Product Name	Storage	Cat. No.	PKG size	Price (US\$)
Bullet Blocking One for Western Blotting	4°C	13779-56	50 ML	\$25.00
	4°C	13779-14	200 ML (50 ML x 4)	\$60.00
	4°C	13779-01	1 L	\$180.00